

There's the Professional Solution. Then There's Everything Else.



Medical professionals need professional tools. For point-of-care glucose testing, leave the at-home strip meters on the shelf. Choose the professional-grade HemoCue® Glucose 201 System.

 HEMOCUE®

HemoCue® Glucose 201 System

The professional solution for lab-quality glucose testing at the point-of-care.




Ideal for both patients with diabetes and those at-risk, the HemoCue® Glucose 201 System provides lab-quality glucose from a fingerstick, or venous sample.

Unlike glucose strip meters, the HemoCue® Glucose 201 System can be used to support the clinical evidence in the diagnosis and treatment of diabetic patients. Providing the ability to perform Oral Glucose Tolerance Tests (OGTT), the HemoCue® Glucose 201 System is the professional solution for in-office testing.

Early Intervention Means Everything

Diabetes is an epidemic in the United States. Early and accurate diagnosis is vital to giving patients a fighting chance against diabetes.¹ Why trust at-home strip meter tests for something so important?

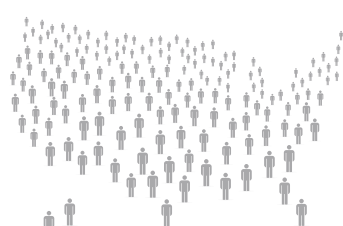

34.2
MILLION
PEOPLE HAVE
DIABETES

1  **5**
OUT OF
DON'T KNOW
THEY HAVE IT

EVERY
8 **SECONDS**
SOMEONE
GLOBALLY
DIES FROM
DIABETES



Sources:


<https://www.cdc.gov/diabetes/pdfs/data/statistics/national-diabetes-statistics-report.pdf>

<https://www.cdc.gov/diabetes/basics/quick-facts.html>

IDF Diabetes Atlas, 9th Edition (2019), pg. 54; https://www.diabetesatlas.org/upload/resources/2019/IDF_Atlas_9th_Edition_2019.pdf

Professional-grade Accuracy Made Easy

The HemoCue® Glucose 201 System includes the analyzer and microcuvettes. Designed for professional use only, our system delivers the results you can rely on to provide the best care possible.

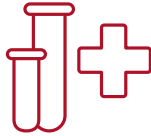

Lab-quality results cleared for aiding in the diagnosis of diabetes; provides the ability to perform CLIA-waived OGTT in-office²




Can be used on patients at-risk for diabetes³



Remain compliant with Clinical Laboratory Improvement Amendments (CLIA) regulations

HemoCue OnCue™ Support

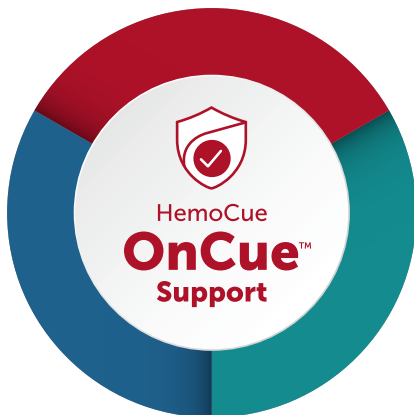

HemoCue OnCue™ Support is a comprehensive program designed to meet your needs as quickly and conveniently as possible.

Powered by HemoCue employees, you get answers and education from the experts who know our products best.


At your office, online, or on the phone, we're ready to do more for you.

¹Brunton SA, Davis SN, Renda SM. Early intervention to achieve optimal outcomes in type 2 diabetes: a case presentation. Clin Cornerstone. 2006;8 Suppl 2:S6 S18. Review. PubMed PMID: 16939876.


²Glucose product insert ver. 150702 180514

³https://care.diabetesjournals.org/content/42/Supplement_1/S13


Get Results in 3 Steps

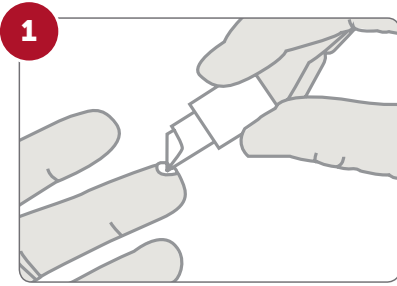

1 Fill the microcuvette.

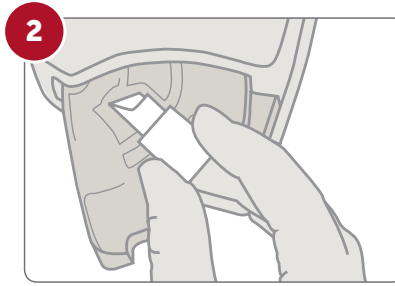

2 Place the microcuvette into the analyzer.

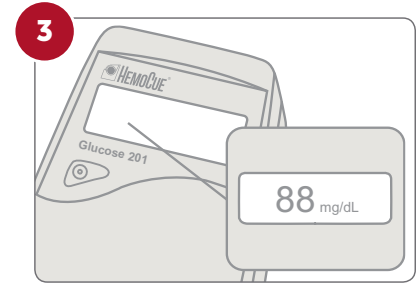

3 View results.

HemoCue® Glucose 201 System Product Specifications

Principle:	Modified glucose dehydrogenase in which the total amount of glucose is measured at the end point photometrically
Calibration:	Factory calibrated and traceable to ID GC-MS method, needs no further calibration and no coding
Sample Material:	Capillary, venous or arterial whole blood including neonatal blood
Displayed Range:	Plasma equivalent values 0-444 mg/dL
Results:	40-240 seconds
Sample Volume:	Approximately 5 µL whole blood
Shelf Life:	Vials and each individually wrapped microcuvette: 9 months from manufacturing
Storage Temperature:	Analyzer: 32-122 °F (0-50 °C); Microcuvette: below 46 °F (8 °C), incl. freezer. Up to 3 days in room temperature 59-86 °F (15-30 °C)
Operating Temperature:	64-86 °F (18-30 °C)
Operating Humidity:	Up to 90 % relative humidity, non-condensing
Power:	AC adapter or 4 AA batteries
Interface:	Connectable to HemoCue® DMU Docking Station
Quality Control:	Built-in self-test; system can be verified using liquid control
CPT Code*:	82947-QW: Blood glucose (sugar) level 82950-QW: Blood glucose (sugar) level after receiving dose of glucose 82951-QW: Blood glucose (sugar) tolerance test, 3 specimens 82952-QW: Blood glucose (sugar) tolerance test, each additional beyond 3 specimens

Product Order Number

120706

HemoCue® Glucose 201 System

**Choose the professional solution – the HemoCue® Glucose 201 System.
Contact your HemoCue or preferred distribution representative today.**

*The CPT codes provided are based on AMA guidelines and are for informational purposes only. CPT coding is the sole responsibility of the billing party. Please direct any questions regarding coding to the payer being billed.



HemoCue America | California, USA | Phone 800.881.1611 | Fax 800.333.7043
web@hemocue.com hemocue.us
©2021 HemoCue America. All rights reserved. MMUS-0435[3]